

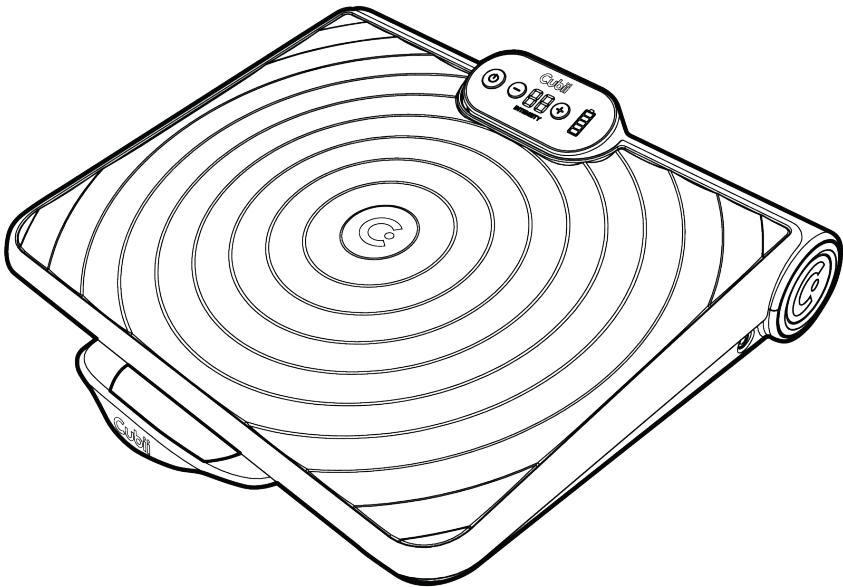


Cubii® **Revive™**

SEATED VIBRATION FOOTPLATE

Owner's Manual

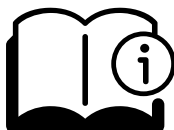

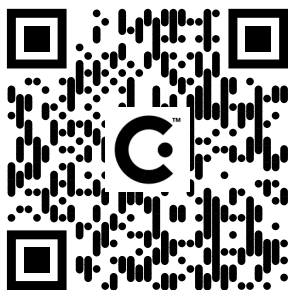

Additional product information and manuals can be found online

Des renseignements supplémentaires sur les produits
sont disponibles en ligne

Puede encontrar información adicional sobre el producto en línea

manuals.cubii.com

TABLE OF CONTENTS

EN

Before You Begin..... 2

Safety Warnings & Instructions

Consult Your Physician

Getting Started..... 5

Product Size and Weight

What's in the Box

Warning Labels

Serial Number Location

Using Your Cubii Revive..... 6

Positioning the Unit

Operating & Adjusting the Intensity

Charging & Checking the Battery

General Information..... 8

Return Policy & Warning

Troubleshooting

Care and Maintenance

Important Health & Safety Instructions

Regulatory Notices

Welcome to the Cubii Communitii®!

Congratulations on taking the next step on your wellness journey and bringing the Cubii Revive into your routine. The Cubii Revive will help rejuvenate the muscles in your lower body, and we are excited for you to get started! Before you begin, we want you to be familiar with important information on how to use the Cubii Revive safely so that you and your loved ones get the most out of your experience.

Let's review the safety notes and warnings now.

If you have any questions, please contact:

+1-844-282-4401

support@mycubii.com


cubii.com/help

BEFORE YOU BEGIN

Safety Warnings & Instructions

WARNING Consult your physician before starting this or any exercise program. This is especially important if you have never exercised before, are pregnant, have a pacemaker, or suffer from any health problems. If at any time during exercise you feel faint, dizzy, or experience pain, stop and consult your physician.

WARNING Do not use when standing. Only use when sitting. Failure to follow all warnings and instructions could result in serious injury or death.

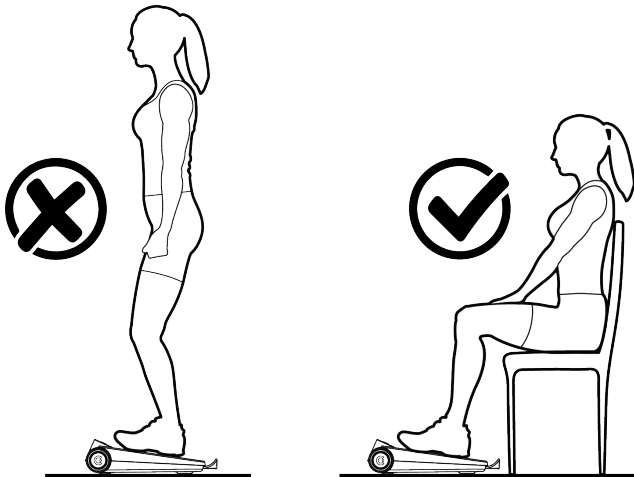

WARNING Use caution when standing after use. For best results, allow a short rest before standing.

WARNING To reduce the risk of serious injury, read the following Safety Instructions before using Cubii Revive.

Read all warnings and cautions posted on the Cubii Revive unit.



1. Cubii Revive should only be used after a thorough review of the Owner's Manual.
2. Do not subject Cubii Revive to heavy shock or treat it with excessive force.
3. If the unit is damaged or not working properly, contact Fitness Cubed, Inc. for service or repair. DO NOT open the unit to service it yourself.
4. To reduce the risk of burns, fire, electric shock or injury to persons, the unit should never be left unattended when plugged in. Unplug the unit when not in use.
5. The Cubii Revive is not waterproof.
6. Do not use cleaners or solvents to clean Cubii Revive.
7. Do not use outdoors or expose to extreme high or low temperatures. Do not leave Cubii Revive in direct sunlight for an extended period.
8. Do not operate under a blanket or pillow. Excessive heating can occur and may cause fire, electric shock or injury.
9. Do not operate around aerosol products or when oxygen is being administered.

10. Close supervision is necessary when the unit is used by children or people with disabilities.
11. Cubii Revive should be used by only one person at a time.
12. Do not carry the unit by its cord or use the cord as a handle.

Consult Your Physician

It is always important to consult your physician before starting an exercise program. This is particularly true if any of the following apply to your current medical condition:

- Chest pain or pain in the neck and/or arm
- Shortness of breath
- A diagnosed heart condition
- Joint and/or bone problems
- Currently taking cardiac and/or blood pressure medications
- Have not previously been physically active
- Dizziness or blurred vision
- Those with physical handicaps

If you feel any of the physical symptoms listed above when you start your exercise program, contact your physician right away.

GETTING STARTED

Product Size and Weight

The assembled Revive measures:



16.9" L x 14.6" W x 5.1" H (43cm x 37cm x 13cm)

The weight of the unit is approximately 8.8lb (4kg).

What's in the Box

Revive Platform


Power Supply Cord

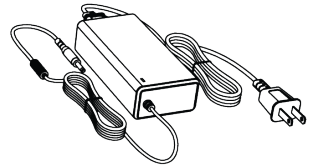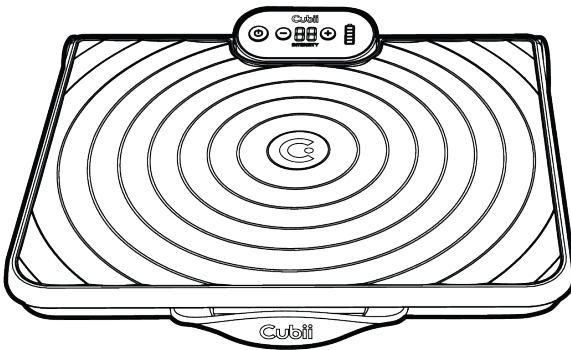

Serial Number Location

The serial number is located on your Cubii Revive product label and the content panel of the box.

Record your serial number here for future reference:

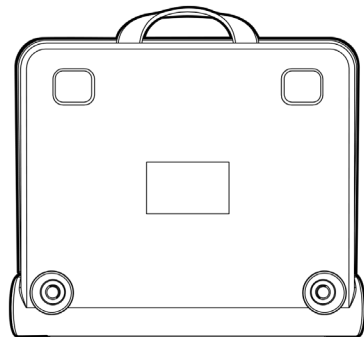


USING YOUR CUBII REVIVE

Positioning

The Cubii Revive can be used either in an upright position, such as with an office chair, or in a more reclined position, such as when used on a sofa.


- Cubii Revive occupies a floor area of 16.9" w x 14.6" L (43cm x 37cm). Please leave a clear zone around the training area that is free of any obstructions.
- Use Cubii Revive indoors on a level surface. Placing the unit on carpeting or a rug is recommended. Keep it away from moisture and dust.
- Do not place Cubii Revive in a location where it may cause an obstruction.

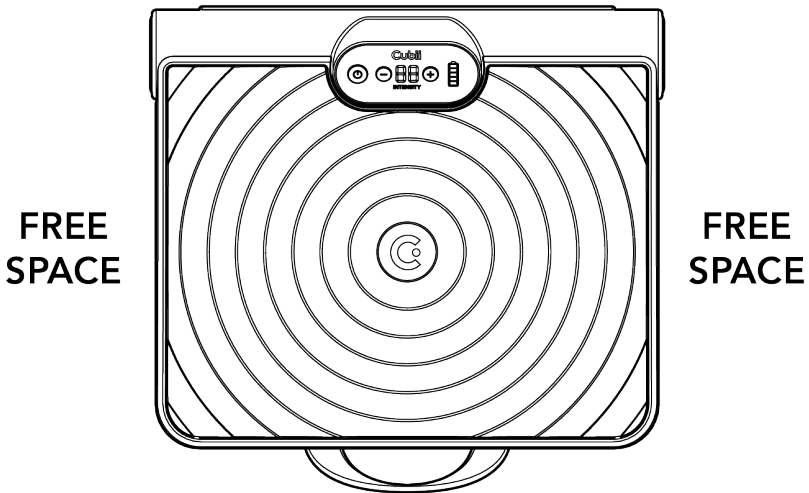

Operating and Adjusting the Intensity

To operate the Cubii Revive, follow the steps below:

1. Place the Cubii Revive on a stable surface in front of a chair or sofa. Keep both feet firmly and securely on the platform while using the Cubii Revive.
2. Press the power button to turn on your Cubii Revive, and then place your feet (socked feet or with shoes on) onto the foot pad.
3. Adjust the intensity between Levels 1 - 10 using the "+" or "-" buttons on the display.
4. To stop using the Revive, press the power button to turn the unit off.

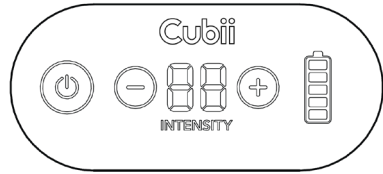

Recommended Usage: We recommend beginning at Level 1 and increasing the intensity gradually over time. Please note that the closer feet are to the center of the foot pad, the more intense the vibrations will feel. We recommend using the Cubii Revive up to 15 minutes at a time. It can be used several times a day

Charging and Checking Battery

The Cubii Revive comes equipped with a rechargeable lithium-ion battery. The unit can operate with or without being plugged in. Charging your Cubii Revive completely can take up to 7 hours. If the unit needs to be charged, follow the steps below:

- Before plugging in the unit, check the power supply cord and be sure that every part is in good working condition. All parts must be firmly assembled.
- Plug the unit in to charge. Plug is located on the right side of the unit.
- When the unit is charging, there should be no obstruction on the cord (no pitching, pulling, etc.).
- Check the battery level at any time on the built-in-screen. The number of bars in the battery symbol indicates the charge left in the battery. When the battery has 15 minutes of charge left, the bottom light will blink. The battery indicator will remain solid when charging is complete but still plugged in.

GENERAL INFORMATION

Return Policy and Warranty

For more information about how to register your unit and our warranty and return policy, scan the QR code at the beginning of the manual, call our Customer Service number +1-844-282-4401 or visit Cubii.com.

Troubleshooting

For troubleshooting and assistance with your Cubii Revive, you can email us at: support@mycubii.com you can visit our help center by visiting: cubii.com/help

Care and Maintenance

Special attention should be given to the following:

Periodically wipe the unit down with a damp cloth, paying special attention to the rubber padding on the platform. Only Fitness Cubed, Inc. supplied components shall be used to maintain/repair Cubii Revive.

Regulatory Notices

1. This device complies with Part 15 of FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) This device must accept any interference received, including interference that may cause undesired operation.

2. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

RF warning statement:

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure conditions without restriction.

Li-ion Rechargeable Battery

Rated: 11.1V 2200mAh 24.42Wh

SHENZHEN MOTTCELL NEW ENERGY TECHNOLOGY CO., LTD

Warning: Risk of fire and burns. Do not open, crush, heat above 45°C or incinerate.

Follow manufacturer's instructions.

® TRADEMARK OF FITNESS CUBED INC. REGISTERED IN THE US AND CANADA.

CAN ICES-003 (B) / NMB-003 (B)

This device contains license-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's license-exempt RSS(s).

Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

L'émetteur/récepteur exempt de licence contenu dans le présent appareil est conforme aux CNR d'Innovation, Sciences et Développement économique Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes :

- 1) L'appareil ne doit pas produire de brouillage;
- 2) L'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement

Model No: VB01

Net Weight: 9.5lb (4.3kg)

Gross Weight 13.4lb (6.1kg)

Product Size: 16.9"L x 14.6"W x 5.1"H (43cm x 37cm x 13cm)

Packaged Size: 19.3"L x 17.1"W x 6.6"H (49.2cm x 43.5cm x 16.9cm)

Maximum Bearing Weight: 300lbs. (136.08kg)

Durability Class: HC

Designed in Chicago, IL

320 W. Ohio St. Suite 700, Chicago, IL 60654

Made in China





Customer Support
support@mycubii.com
cubii.com/help
+1-844-282-4401